



# Real-time image denoising of mixed Poisson–Gaussian noise in fluorescence microscopy images using ImageJ

VARUN MANNAM,<sup>1,\*</sup>  YIDE ZHANG,<sup>1,2</sup>  YINHAO ZHU,<sup>1</sup> EVAN NICHOLS,<sup>3</sup> QINGFEI WANG,<sup>3</sup> VIGNESH SUNDARESAN,<sup>4</sup> SIYUAN ZHANG,<sup>3</sup> CODY SMITH,<sup>3</sup>  PAUL W. BOHN,<sup>4,5</sup> AND SCOTT S. HOWARD<sup>1,6</sup> 

<sup>1</sup>Department of Electrical Engineering, University of Notre Dame, Notre Dame, Indiana 46556, USA

<sup>2</sup>Caltech Optical Imaging Laboratory, Andrew and Peggy Cherng Department of Medical Engineering, California Institute of Technology, Pasadena, California 91125, USA

<sup>3</sup>Department of Biological Sciences, University of Notre Dame, Notre Dame, Indiana 46556, USA

<sup>4</sup>Department of Chemical and Biomolecular Engineering, University of Notre Dame, Notre Dame, Indiana 46556, USA

<sup>5</sup>Department of Chemistry and Biochemistry, University of Notre Dame, Notre Dame, Indiana 46556, USA

<sup>6</sup>e-mail: showard@nd.edu

\*Corresponding author: vmannam@nd.edu

Received 10 November 2021; revised 1 February 2022; accepted 6 February 2022; published 22 March 2022

Fluorescence microscopy imaging speed is fundamentally limited by the measurement signal-to-noise ratio (SNR). To improve image SNR for a given image acquisition rate, computational denoising techniques can be used to suppress noise. However, common techniques to estimate a denoised image from a single frame either are computationally expensive or rely on simple noise statistical models. These models assume Poisson or Gaussian noise statistics, which are not appropriate for many fluorescence microscopy applications that contain quantum shot noise and electronic Johnson–Nyquist noise, therefore a mixture of Poisson and Gaussian noise. In this paper, we show convolutional neural networks (CNNs) trained on mixed Poisson and Gaussian noise images to overcome the limitations of existing image denoising methods. The trained CNN is presented as an open-source ImageJ plugin that performs real-time image denoising (within tens of milliseconds) with superior performance (SNR improvement) compared to conventional fluorescence microscopy denoising methods. The method is validated on external datasets with out-of-distribution noise, contrast, structure, and imaging modalities from the training data and consistently achieves high-performance (>8 dB) denoising in less time than other fluorescence microscopy denoising methods. © 2022 Optica Publishing Group under the terms of the


Optica Open Access Publishing Agreement

<https://doi.org/10.1364/OPTICA.448287>

## 1. INTRODUCTION

In modern biology, fluorescence microscopy plays a vital role in functional and structural imaging [1]. However, the imaging speed of fluorescence microscopy is limited by the inherent noise in the system [2]. Fluctuations in the photon emission events contribute to quantum noise/shot noise (from the stream of photons), which follows Poisson statistics per time period. Additionally, microscope electronic components contribute Johnson–Nyquist noise (or thermal noise), which follows a Gaussian distribution per time period. The combination of both of these contributions leads to fluorescence microscopy systems containing a mixed Poisson–Gaussian (MPG) noise. Due to this inherent noise in the system, the acquired image has a limited signal-to-noise ratio (SNR) that hinders the underlined ground truth information about the biological sample. Therefore, a robust denoising approach is required to address these two distinct noise distributions.

In fluorescence microscopy, the statistics of each noise source affect images differently. To illustrate these differences, we simulate the effects of Gaussian noise and Poisson noise individually on a clean image in Fig. 1. The clean image represents bovine pulmonary artery endothelial (BPAE) cells (Invitrogen FluoCells slide#1 F36924, mitochondria labeled with MitoTracker Red CMXRos, F-actin labeled with Alexa Fluor 488 phalloidin, nuclei labeled with DAPI) imaged with a commercial fluorescence microscope (Nikon A1R-MP laser scanning confocal microscope equipped with a Nikon Apo LWD 40×, 1.15 NA water immersion objective and more details are provided in Supplement 1 Note S1) with a 4 mW excitation power and 50 frame averaging to produce a high SNR image. Gaussian noise and Poisson noise are added to the clean image (SNR of 25 dB with reference to the clean image) and shown as “With Gaussian noise,” and “With Poisson noise,” images. The “Clean image” is obtained by averaging 50 noisy images in the same field of view (FOV). Here the Gaussian noise

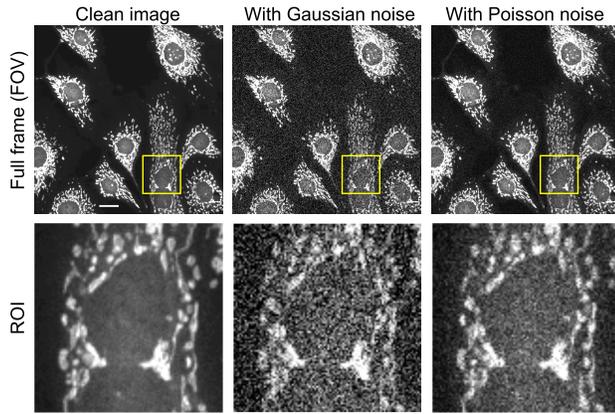

**Fig. 1.** Illustration of different types of noise with a BPAE cell sample. The top row indicates full frame images, and the PSNR of the synthetic noisy image is 25 dB with reference to clean image (obtained by averaging 50 noisy images in the same FOV). The bottom row shows the selected region of interest (ROI) as shown in the yellow box of the corresponding top row images. Full frame and ROI dimensions are  $512 \times 512$  and  $125 \times 125$ , respectively. Scale bar: 20  $\mu\text{m}$ .

is independent of clean image pixel intensity, whereas the Poisson noise is pixel intensity dependent. We can observe this behavior in the marked region of interest (ROI) center region, i.e., the smooth gray area (nuclei region). Hence, Poisson noise and Gaussian noise affect the clean image differently, which makes MPG noise denoising a challenging task.

To improve fluorescence microscopy SNR, one can increase the sample dosage by increasing either the excitation power or pixel dwell time. However, this method could lead to photobleaching or photodamage and cannot be used for imaging real-time dynamics in living animals. Furthermore, lower dosage leads to noisy (with low SNR) images, thereby posing a challenge for applications such as cell classification and segmentation [3]. Denoising algorithms are available to improve image SNR, but either are designed specifically for Gaussian or Poisson noise (thereby not suitable for other distribution or MPG noise) or are capable of denoising MPG noise but require long computation times to achieve relatively small SNR improvement. Therefore, there is a need for a computationally fast single-frame denoising method that works for all experimental noise distributions: Poisson, Gaussian, and MPG.

Reiterating the mathematical notation of different noise distributions using first and second order moments (or simply mean and variance) [4], the Gaussian noise distribution is represented as standard normal distribution  $n_g \sim \mathcal{N}(0, b)$ , where the mean is zero and the variance is  $b$ . Similarly, the Poisson noise distribution can be represented as  $n_p \sim \mathcal{N}(y_i, ay_i)$ , where the true pixel value ( $y_i$ ) is the mean and Poisson variance is  $ay_i$  (conversion factor or quantum efficiency of  $a$  times  $y_i$ ). Assuming that the Poisson and Gaussian processes are independent, the MPG noise distribution is represented as  $n_{\text{mixed}} \sim \mathcal{N}(y_i, ay_i + b)$ , where the mean is the true pixel value ( $y_i$ ) and the variance is the addition of two independent noise variances ( $ay_i + b$ ). Hence, MPG noise distribution is similar to Poisson noise with an additional Gaussian noise variance ( $b$ ).

In addition to the above-described noise statistics, fluorescence microscopy systems using detector arrays [e.g., charge-coupled device (CCD) or complementary metal-oxide semiconductor (CMOS)] can experience a “fixed-pattern noise” due to differences in detector responsivity and dark current between pixels. Such

noise commonly appears as horizontal or vertical “stripes” in low-signal regions of an image. The signal variance in the presence of fixed-pattern noise also follows the form  $ay_i + b$ , but with  $a$  and  $b$  varying at each pixel. Denoising such systems using calibration or noise statistics is possible [5], but requires prior details on each specific system. Fixed-pattern noise can also be addressed using a photon event centroid algorithm [6], flat-field correction method [7], and automatic correction of CMOS-related noise using camera properties combined with sparse filtering [5], or calibration by measuring the pixel-specific dark current and gain, then compensating for it via subtraction and division, respectively. In this paper, we will primarily consider images that do not experience fixed-pattern noise or have compensated for fixed-pattern noise via calibration.

To increase the SNR of an image, additional information must be included to estimate the original, denoised image. If the noise distribution is known *a priori* to be Gaussian, simple filtering (by assuming noise exceeds signal at high spatial frequencies by using mean, median, bilateral filters) or machine learning (ML)-based approaches such as non-local means (NLM) [8], block matching 3D filtering (BM3D) [9], K-means singular value decomposition (K-SVD) [10], expected patch log-likelihood (EPLL) [11], and weighted nuclear norm minimization (WNNM) [12] can all be used to leverage knowledge of the noise distribution and to leverage information from similar features within the same image. If the noise distribution is known to be Poisson, the above methods can be used by (1) converting the Poisson noise into Gaussian noise using a variance stabilization transform (VST) [e.g., generalized Anscombe transform (GAT [13])], (2) denoising using popular Gaussian denoising algorithms, or (3) converting the denoised results back using an inverse VST [13]. Such popular Poisson denoising methods include Poisson NLM (PNLM) [14], non-local principle component analysis (NL-PCA) [15], VST + NLM and VST + BM3D, which apply VST, NLM or BM3D Gaussian denoising methods, and inverse VST, sequentially, on an image. If the noise is from a MPG distribution, more complicated methods are needed to estimate the relative contributions of each process at each pixel (e.g., Poisson unbiased risk estimate – linear expansion of thresholds (PURE-LET) [16]).

Overall, all of these image denoising methods have trade-offs among denoising performance, computation speed, and applicability to a specific noise distribution. In this paper, we show how a convolutional neural network (CNN) trained on MPG noise can in fact achieve high-performance denoising, quickly, on images with Poisson and/or Gaussian noise distributions.

Recently, CNN-based fluorescence microscopy image denoising has demonstrated improved denoising performance and speed [17]. CNN-based image denoising methods can be classified into two categories: supervised and self-supervised. Supervised methods require training the network on a dataset prior to performing denoising on experimental images and include content-aware image restoration (CARE) [17], denoising CNNs (DnCNN) [18], and Noise2Noise [19]. Self-supervised methods train and perform denoising using the same experimental image and include Noise2Void (N2V) [20], probabilistic N2V (PN2V) [21], Noise2Fast [22], and other methods [23,24]. Supervised CNN training methods typically provide higher performance than self-supervised methods due to the inclusion of additional information in target images. In some conditions, self-supervised denoising methods can yield performance very comparable to

supervised methods [19]. In this paper, we utilize a supervised CNN model that is capable of removing MPG noise from a single frame nearly in real time (within a few ms). By producing and utilizing a large, multiple species, and multiple microscope modality training dataset [25], the trained model outperforms existing image denoising methods. More importantly, our method demonstrates superior denoising performance compared to the existing image denoising approaches including ML-based methods independent of noise distribution (either Poisson or Gaussian or MPG noise) presented in most fluorescence microscopy images.

To address the fundamental problem of denoising microscopy images containing MPG noise sources, we train and evaluate well-known CNN models on MPG noise microscopy data from the Fluorescence Microscopy Denoising (FMD) dataset [26] in Sections 2 and 3. In Section 4, the trained ML model is validated by evaluating performance in a diverse set of microscopy images, including the Widefield2SIM(W2S) dataset [27], with out-of-training-distribution SNR and contrast. The method consistently produces high-quality denoised images and is packaged as an ImageJ plugin to reduce the barrier to entry for applying denoising CNNs to microscopy data.

## 2. METHODS

As described in Section 1, the fundamental stochastic nature of MPG noise requires introducing some *a priori* knowledge to the denoising process. Analytical approaches use knowledge of the statistical properties of noise, but are computationally slow and not well suited for MPG sources. Deep CNNs, however, are data driven and expressive enough to overcome the complex statistical nature of MPG noise. To perform denoising using CNNs, a network architecture and training datasets must be chosen. In this paper, we choose to focus on the Noise2Noise [19] and DnCNN network [18] architectures for two reasons. First, they have demonstrated superior performance to other CNN architectures [20]. Second, robust performance in diverse microscopy environments requires training on large datasets. Supervised training of the Noise2Noise and DnCNN models can be achieved in reasonable time ( $\approx 4$  h for the entire FMD dataset on a single Nvidia 1080-ti GPU) while training other networks (e.g., N2V) on the dataset can take a prohibitively long time ( $\approx 3.50$  h for a single image training on the GPU). N2V is a popular denoising method that is capable of self-supervised learning, so we will therefore compare supervised trained model performance to self-supervised N2V performance in Section 4.A. We chose to train these two models on our FMD dataset [25,26] that contains 12,000 raw noisy images from confocal [28], two-photon [29], and widefield [30] microscopes. The FMD dataset also helps to benchmark various denoising techniques on the fluorescence images. The creation of the training and testing datasets from the FMD dataset is explained in Supplement 1 Note S1.

The first ML model we evaluate is based on a deep CNN with the Noise2Noise architecture [19] trained to denoise fluorescence images with MPG noise. The ML model is similar to the originally published Noise2Noise architecture [19] with the addition of a nonlinear activation (tanh) layer after the final convolutional layer [26], adaptive learning (one cycle policy) during the training phase, and carefully tuned hyper-parameters such as batch size and learning rate for achieving best fluorescence microscopy image denoising performance. Here onwards, we call the modified model “Noise2Noise plugin.”

In the Noise2Noise plugin CNN model, the input and target are both noisy images within the same FOV (which is the case when the ground truth is absent or its extraction is difficult). This ML model is beneficial compared to our second ML model, which is based on the DnCNN architecture [18] that requires ground truth since getting the clean image (typically acquired with the increased dosage) is usually difficult, especially when imaging living species. Figure 2 shows the complete Noise2Noise plugin architecture as our first ML model. Here we divide the ML architecture into five blocks. Blocks A to C represent an encoder structure, and blocks C to E represent a decoder structure. In the encoder, the image feature size ( $256 \times 256$ ) is reduced by selecting only essential features of the image (with more channels to select different features), and at block C, the images with feature map size  $8 \times 8$  are less noisy when compared to the noisy input image. In the decoder, this noise-free image is restored to the original noisy image feature size ( $256 \times 256$ ). In this work, both input and target are noisy images. At the  $i$ th pixel of the input and target images, the noisy measurements are  $z_1 = y + n_1$  and  $z_2 = y + n_2$ , where  $y$  is the real value and  $n_1$  and  $n_2$  are MPG noise values corresponding to the input and target, respectively. The Noise2Noise plugin model estimates the real value  $y$  instead of the noisy target  $z_2$  since  $n_2$  is random. The Noise2Noise plugin model can predict the clean image since training happens on the noisy image pair as input and target, and this ML model works well as long as the clean image is the average of noisy target images within a FOV. The other ML model is based on a deep CNN with the DnCNN architecture trained to denoise the fluorescence images with MPG noise. The DnCNN ML model can estimate the residual (only noise) for a given noisy input image, and the denoised image is obtained by subtracting the estimated residual from the noisy input image. Unlike the supervised DnCNN ML model that requires clean images as the target, the Noise2Noise plugin ML model utilizes the noisy image as the target. The training and testing of the proposed ML models are explained in Supplement 1 Note S2.

Interestingly, both the Noise2Noise and DnCNN in literature use relatively simple U-Net and fully CNNs, respectively, yet achieve high-quality denoising results. However, more advanced architectures and methods are commonly required to perform microscopy image identification and segmentation [31–34]. While the model capacity of simple U-Net and CNNs is reduced compared to advanced methods, denoising models require training a model only to discern between the presence and absence of general common features, while image identification and segmentation models must be trained to discriminate between specific high-detail features. We therefore employ simple, yet high-performance, denoising networks in this study and evaluate performance on a variety of image samples in Section 4.

While the proposed ML approach is fast and accurate, the complexity of implementing the ML model is a barrier to entry. To address this barrier, we integrated our ML models as plugins to ImageJ, one of the most widely used biomedical image processing packages [35]. Two ImageJ plugins are presented for two ML models: Noise2Noise plugin and DnCNN plugin. We have developed our ImageJ plugins based on the ImageJ-tensorflow library [36]. Here the library provides the minimal set of functions in ImageJ to load the pre-trained ML model weights as a graph and perform a set of convolutional layer operations on the input image (loaded by ImageJ), which results in the output of the trained ML model. To accommodate more image formats, we




due to an offset in the noisy image. More details can be found in [40] Section II (A). For these cases, we set the estimated Gaussian variance to zero. From Table 1, we also observe that the estimated Poisson noise variance ( $a$ ) in confocal and two-photon microscopy dominates the Gaussian noise variance ( $a \gg b$ ). In the case of widefield microscopy, both noise parameters exist ( $a, b > 0$ ), which indicates the mixture of Poisson–Gaussian noise. Also, the Gaussian variance ( $b$ ) is in the range of  $\sim 2e-4$ , whereas the Poisson variance ( $a$ ) varies depending on the signal intensity. This behavior is observed in the widefield microscope, where the  $a$  value of BPAE (F-actin) is much larger when compared to the  $a$  value of BPAE (mitochondria) and BPAE (nuclei). The estimation of these noise parameters, i.e., the Poisson variance ( $a$ ) and the Gaussian variance ( $b$ ), shows a clear map between different noises in each microscope modality. Hence, we broadly divide the microscopy modalities into two categories: first, confocal and two-photon microscopy with Poisson-dominated noise; second, widefield microscopy with MPG noise. For each class of noise type, we show the trained model performance on a single noisy image and the statistics on the complete test-mix dataset.

In this work, we first select a single test image with Poisson-dominated noise, which is captured by a confocal microscope, and another test image with MPG noise, which is captured by a widefield microscope. Here these images are considered as specific images that do not reflect the complete test dataset.


Figure 4 illustrates a zebrafish embryo image with Poisson-dominated noise [EGFP labeled Tg (sox10:megfp) zebrafish at two days post fertilization] captured with a commercial confocal microscope (more details are provided in Supplement 1 Note S1). Similarly, Fig. 5 shows an image of a BPAE cell (slide #1, F36924 containing nuclei, F-actin, and mitochondria) with MPG noise captured by a commercial widefield microscope. In these figures, the target intensity images were generated by taking the average of 50 noisy images in the same FOV. We observe that our ML-based Noise2Noise plugin denoising method reduces both Poisson-dominated and MPG noise significantly by providing a PSNR improvement of 10.16 dB and 10.5 dB in zebrafish and BPAE cells, respectively.

To quantitatively evaluate whether the denoised images contain similar image features as the clean image, we calculate the structural similarity index measure (SSIM) between the two. The SSIM compares luminance, brightness, and contrast values as a function of position [41]. SSIM measures the similarity between two images on a scale of zero to one, with one being perfect fidelity. In addition, SSIM can be correlated to the PSNR, and more details are provided in [42]. From the images in Fig. 5, the denoised images appear visually similar to the respective target images, and with denoised images, SSIM values of zebrafish and BPAE samples are 0.84 and 0.83, respectively.

Second, the PSNR distribution of the raw images from the test-mix data and the corresponding denoising results using our ML models are provided in Fig. 6. Figure 6 shows the trained Noise2Noise plugin ML model (with an ILR of  $5e-4$ ) image denoising performance with complete statistics on the test-mix dataset raw images (48) from our FMD dataset. Also, for image denoising, we train another ML model (requires a clean target), which is a DnCNN ML model that can estimate the residual (only noise) for a given noisy input image. Here the denoised image is obtained by subtracting the estimated residual from the noisy input data. Figure 6 also shows complete statistics on the test-mix

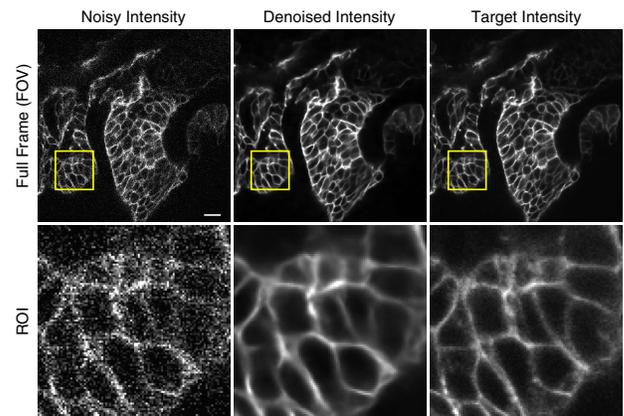


**Fig. 4.** Image denoising results using the Noise2Noise plugin on a fixed zebrafish embryo [EGFP labeled Tg(sox10:megfp) zebrafish at two days post fertilization] captured with confocal microscopy [pixel dwell time of  $2 \mu\text{s}$  and pixel width of 300 nm, with 10% of power and PMT gain of 140 (model number: Hamamatsu, PMT R10699); excitation and emission wavelengths are 488 nm and 509 nm, respectively]. The top row indicates the full-frame (of size  $512 \times 512$ ) of noisy input, denoised output, and target, and the bottom row indicates the region of interest (ROI: marked in the yellow square of size  $100 \times 100$ ) from the respective top row images.

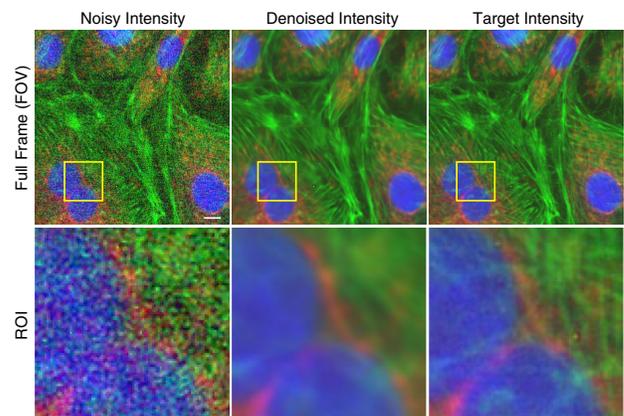

**Fig. 5.** Image denoising results using the Noise2Noise plugin on the BPAE cell [labeled with MitoTracker Red CMXRos (mitochondria), Alexa Fluor 488 phalloidin (F-actin), and DAPI (nuclei) with excitation/emission wavelengths of 561/599, 488/512, and 405/461, respectively; Invitrogen FluoCells prepared slide#1: F36924] captured with widefield microscopy [exposure time of 200 ms, frame rate of 5 Hz, high definition color camera and controller (DS-Fi1 and DS-U2) with CCD gain of 46 and pixel width of 170 nm]. The top row indicates the full-frame (of size  $512 \times 512$ ) of noisy input, denoised output, and target, and the bottom row indicates the region of interest (ROI: marked in the yellow square of size  $100 \times 100$ ) from the respective top row images. We imaged three times for the same FOV, each time with a different filter block (DAPI for nuclei, FITC for F-actin, TRITC for mitochondria), to acquire the multi-channel (color) fluorescence image of the cells.

data using the DnCNN architecture ML model with the same ILR of  $5e-4$ . The ML model with Noise2Noise plugin architecture has a better PSNR when compared to the ML model trained with the DnCNN architecture. Hence, the Noise2Noise plugin model demonstrates superior denoising performance compared to the other ML model, and it does not require the clean target image during the training phase. From the raw data distribution, it is

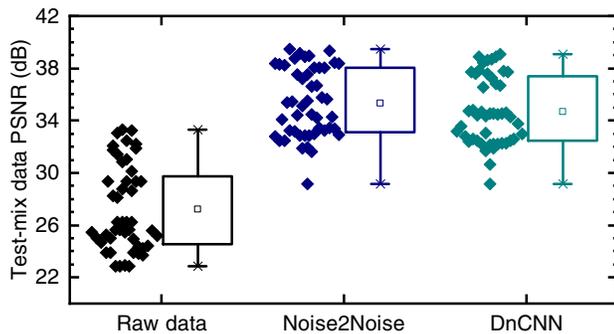


**Fig. 6.** Test-mix dataset (containing 48 images of MPG noise samples from different microscopes in the FMD dataset) raw images PSNR box-plot to compare Noise2Noise plugin and DnCNN plugin ML models' image denoising performances at a fixed initial learning rate (ILR:  $5e-4$ ). Box-plot contains min, max, 25th, and 75th percentile values. Here the square box indicates the mean value. A scattered plot (left to box plot) indicates the true PSNR values used for the statistics.

evident that there exist two distinct groups of data corresponding to Poisson-dominated noise and MPG noise. The average PSNR of the raw images (including Poisson-dominated noise and MPG noise images) and the corresponding denoised images are 27.22 dB and 35.40 dB, respectively. From Fig. 3, the PSNR improvement using the Noise2Noise plugin method is comparable to that of an average PSNR of eight noisy raw images in the same FOV (with PSNR of 36.02 dB), thereby reducing the acquisition time approximately by eight fold.

Finally, we compare the performance of our ML-based image denoising methods with the existing image denoising methods. This PSNR improvement is more substantial than existing image denoising methods such as PURE-LET, VST + NLM, and VST + BM3D, as shown in Fig. 7. Figure 7 shows the existing denoising methods performed on images presented with either Poisson or Gaussian noise. In contrast, our ML-based image denoising methods enable superior denoising performance regardless of the noise distribution or microscopy modality. The PURE-LET, VST + NLM, and VST + BM3D image denoising methods perform well only for the images with Poisson-dominated noise, while the NLM and BM3D image denoising methods perform well only for images with Gaussian noise. In contrast, the ML-based image denoising methods (Noise2Noise plugin and DnCNN plugin) perform well for images with both Poisson-dominated noise and MPG noise.

Our Noise2Noise plugin ML model image denoising method shows significant improvements over existing methods and can provide better results for most fluorescence microscopy image following MPG noise distribution. To increase the impact and integrate it into existing workflows, we package our method as an ImageJ plugin, a popular image analysis tool used by biologists. Our developed ImageJ plugin has two advantages. First, the plugin can perform image denoising independent of microscopy modality since the ML model was trained with MPG noise images with superior performance compared to existing methods. Second, the plugin can denoise in "real time" (e.g., with an image size of  $256 \times 256$  in 80 ms using GPU) for fluorescence microscopy images with MPG noise. Our plugin denoising time also scales accordingly if the noisy image dimensions are different. More details about the ImageJ image denoising plugin are provided in Supplement 1 Note S5. In addition, Supplement 1 Note S3 shows a qualitative and quantitative comparison of the existing image denoising ImageJ plugins (existing and our ML-based image

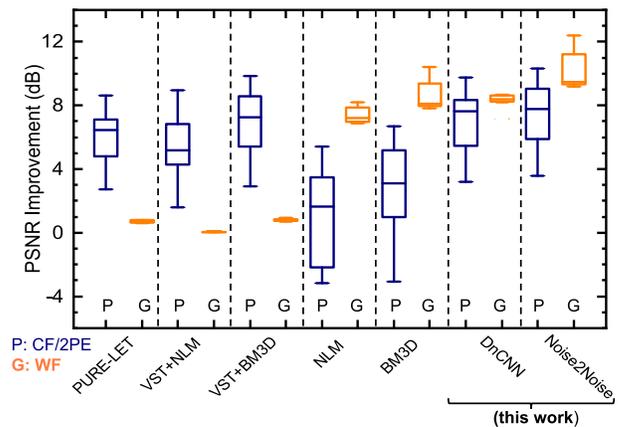

**Fig. 7.** PSNR improvement (denoised image PSNR, raw image PSNR) using different image denoising methods on the raw data with Poisson and MPG noise statistics independently. Here PURE-LET, VST + NLM, VST + BM3D, NLM, and BM3D are existing image methods, and DnCNN plugin and Noise2Noise plugin are our demonstrated supervised trained CNN methods. Blue and orange statistics indicate the box plot results of the PSNR improvement for the Poisson-dominated and MPG noises, respectively. Supervised CNN-based methods outperform compared to existing image methods independent of noise characteristics.

denoising methods) available with the trained Noise2Noise plugin ML model on a test image from the FMD dataset. The image denoising time (for an image of  $256 \times 256$ ) using our Noise2Noise plugin on a CPU (Intel Xeon E5-2680 CPU) and on a GPU (Nvidia GeForce GTX 1080 Ti) are 960 ms and 80 ms, respectively. Therefore, the image denoising speed is faster by 12-fold when using GPU compared to CPU. Overall, the ImageJ plugin can be used to denoise fluorescence images that are obtained with low laser power or at a fast acquisition rate. The source code for image denoising including ML model training and testing in Keras and creation of our ImageJ plugin are provided in [43].

#### 4. MODEL VALIDATION

Model validation outside of the training FMD dataset images and noise distributions is crucial to demonstrate the method's reliability and utility. In the following subsections we validate the Noise2Noise plugin performance on (A) the W2S fluorescence microscopy dataset to compare performance to other analytic and ML-based image denoising methods outside of the FMD training dataset; (B) the W2S dataset to evaluate performance on samples outside of the training dataset the PSNR distribution (out-of-distribution noise) and on a flavin-based auto fluorescent bacteria sample images to evaluate performance outside of the training dataset contrast distribution (out-of-distribution contrast); and (C) the 3D deep residual channel attention network (3D RCAN) dataset to evaluate performance on samples outside of the training dataset structures (out-of-distribution structures). In Section 4.D we list steps to validate the performance on application-specific fluorescence microscopy samples, and in 4.E we discuss the risk and limitations of applying the trained CNN model to higher-dimensional data.

##### A. Denoising Performance on the W2S Dataset

To validate the trained Noise2Noise plugin ML model on images that differ from the FMD dataset, we evaluate denoising plugin

**Table 2. Average PSNR Improvement ( $\Delta$ PSNR = denoised image PSNR, input image PSNR) and Single-Frame Processing Computation Time (Intel Xeon E5-2680 CPU and Nvidia GeForce GTX 1080 Ti GPU) of Noise2Noise Plugin, VST + BM3D, BM3D, and Noise2Void Image Denoising Methods When Applied to 360 Widefield Fluorescence Microscopy Images in the W2S Dataset<sup>a</sup>**


Method	Training	$\Delta$ PSNR (dB)	Time/Frame (s) CPU/GPU
Noise2Noise plugin (this work)	Supervised, FMD	8.25	3.07/0.08
VST + BM3D	Self-supervised	4.37	5.67/None
BM3D	Self-supervised	4.74	4.65/None
Noise2Void (fast settings)	Self-supervised	4.6 (0.7) <sup>a</sup>	15600/400
Noise2Void (default settings)	Self-supervised	4.8 (0.7) <sup>a</sup>	478800/12000

<sup>a</sup>Estimated mean (and standard error of the mean) from 10 random noisy images from the W2S dataset.

performance on 360 widefield fluorescence microscopy images of varying PSNRs from the W2S dataset [27]. We also investigate whether the FMD trained ML plugin model applied to the W2S dataset introduces performance degradation or artifacts compared to three of the best performing self-supervised (no external training) ML methods: N2V [20], VST + BM3D [13], and BM3D [9].

The average PSNR improvement and computation time for the different methods applied to the W2S dataset is given in Table 2. For the Noise2Noise plugin, the trained ImageJ plugin presented here was used; for the VST + BM3D method (in sequential steps: apply VST transform [13], use the BM3D method [9], and apply inverse VST transformation [13]) for image denoising and the BM3D method [9], MATLAB code was used, and for the N2V method, an ImageJ plugin was used [44] with two sets of training parameters to explore performance versus computation time trade-offs (“fast” with 10 epochs and 200 steps, “default” with 300 epochs and 200 steps). In summary, the supervised Noise2Noise plugin trained on the FMD dataset yields higher PSNR improvements compared to self-supervised ML methods. Additionally, the Noise2Noise plugin method requires only convolutions to perform denoising inference, while self-supervised requires re-training (or self-training) on every image. The resulting denoising computation time is therefore significantly shorter for the pre-trained Noise2Noise plugin when computed using either a CPU (Intel Xeon E5-2680 CPU) or GPU (Nvidia GeForce GTX 1080 Ti). While the Noise2Noise plugin, VST + BM3D, and BM3D image denoising methods were performed on the entire dataset to find the mean PSNR improvement over the entire W2S dataset, the N2V method requires a long training time per image, and therefore estimates of the mean improvement (and standard error of the estimate) are presented.

A qualitative comparison of denoising performance is given in Fig. 8 for a representative image from the W2S dataset. Here the target image is an average of 400 noisy images in the same FOV and is provided by the W2S dataset. Additional examples are available in Supplement 1 Note S4 for single-channel and multi-channel test images from the W2S dataset. From these examples, the images qualitatively exhibit lower noise and faithfully represent the target image.

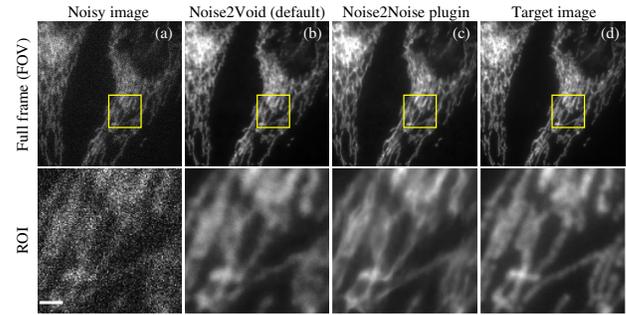

**Fig. 8.** Comparison of the ML-based image denoising plugins in ImageJ. (a) Noisy image, (b) denoised by Noise2Void ImageJ plugin (default training), (c) denoised by the Noise2Noise plugin, and (d) ground truth image of a test image from the W2S dataset. The top row indicates the full frame (of size  $512 \times 512$ ), and the bottom row indicates the region of interest (ROI: marked in the yellow square of size  $120 \times 120$ ) from the respective top row images. Scale bar:  $10 \mu\text{m}$ .


## B. Out-of-Distribution Noise and Contrast

ML models are optimized for inference data distributions that are similar to the training data distribution. The FMD training dataset comprises images with PSNRs in the range of 22–40 dB with contrast ratios [defined as  $(I_{\max} - I_{\min}) / (I_{\max} + I_{\min})$ ] of approximately 0.94. To evaluate model performance outside of the training dataset distribution, we validate the trained Noise2Noise plugin on images of varying PSNRs and contrasts in this section.

The FMD training dataset includes MPG noise from the FMD dataset with PSNRs in the range of 22–40 dB; therefore, the model is expected to perform best when input images with similar noise levels and input images with lower PSNRs outside of the training dataset range can expect lower levels of improvement or artifacts [45]. The trained Noise2Noise plugin model presented here, however, is designed to properly denoise images with out-of-distribution noise using “bias-free” (i.e., no bias term) in the convolutional layers [45]. By not including a bias term in the convolutional layers or batch-normal regularization layers, the Noise2Noise plugin provides PSNR improvements when applied to low-SNR images. In this section, we evaluate denoising performance on the W2S dataset for images with out-of-distribution PSNR.

To validate that the bias-free Noise2Noise plugin model can appropriately denoise low-PSNR out-of-distribution images, we evaluated denoising performance on such data using the method described in [45]: Gaussian noise is synthetically added to images from the W2S dataset to yield images with PSNRs in the range of 5–40 dB, and denoising performance can thereby be evaluated as a function of input PSNR. Results are shown in Fig. 9 (dark blue symbols and line), with significant PSNR improvements ( $>9.6$  dB) for the W2S test data outside of the trained PSNR distribution ( $<22$  dB). For images with the trained noise distribution (input PSNRs of 22–40 dB, light blue shading), the average PSNR improvement is  $\approx 7.06$  dB. PSNR improvement decreases as the input PSNR value increases within the trained noise distribution, and this behavior is consistent with the results shown in Fig. 3.


To validate denoising plugin performance for images with out-of-distribution contrast, we evaluate denoising performance for low-contrast images [for illustration, defined as  $(I_{\max} - I_{\min}) / (I_{\max} + I_{\min})$ ]. This metric defines contrast in the range of zero to one, where one is ideal high contrast. The FMD

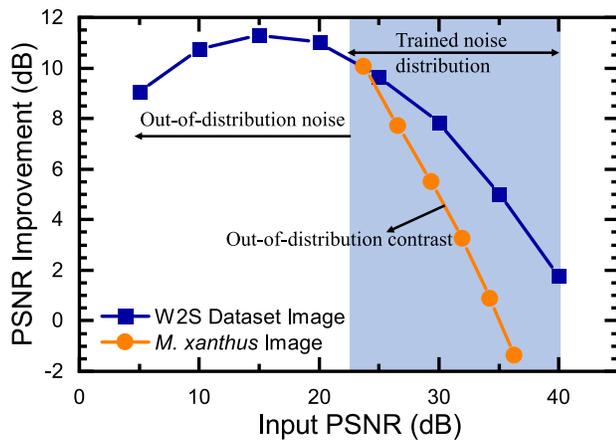


**Fig. 9.** PSNR improvement using the Noise2Noise plugin on test data with acquired PSNR (blue) and contrast (orange) that are out of distribution relative to the PSNR and contrast of the FMD training dataset. Out-of-distribution low-PSNR images were obtained by adding synthetic noise to the W2S dataset. Out-of-distribution contrast images are taken from the *M. xanthus* autofluorescence dataset and averaged (one, two, four, eight, 16, and 32 frames with the same FOV) to obtain images with varying PSNR.

dataset and W2S datasets are both high-contrast datasets with average contrast measures of 0.94 and 0.95, respectively. To evaluate denoising performance for low-contrast fluorescence microscopy images, we prepared an experimental dataset of 600 images (time lapse data) of *Myxococcus xanthus* bacteria auto-fluorescence using widefield microscopy. This sample and imaging method were chosen to produce experimental results with low average contrast (0.52), far outside of the training distribution. Sample preparation and microscope configuration are described in [46,47]: unlabeled bacteria are illuminated with a 458 nm laser to excite flavins in the bacteria, resulting in autofluorescence. The microscope is configured as widefield epi-fluorescence microscopy using a single dichroic mirror as both excitation and emission filters, resulting in the excitation light leakage and/or backscattering being detected by the imaging EM-CCD. The resulting contrast is therefore relatively low due to the low emission intensity from unlabeled autofluorescence and the weak background suppression from using a single excitation blocking filter. A representative result is shown in Figs. 10(a)–10(c), depicting the out-of-distribution noisy *M. xanthus* image, denoised using the Noise2Noise plugin ML model, and target images, respectively. The target image is defined as the average of 600 noisy images in the same FOV. In this example, applying the Noise2Noise plugin on a single noisy frame from this dataset increases PSNR from 23.72 dB to 33.78 dB.

Denoising plugin performance for these low-contrast images can be evaluated as a function of input image PSNR by averaging sequential frames prior to running the denoising plugin. The resulting denoising improvement is given in Fig. 9 (orange symbols). The Noise2Noise plugin achieves high-performance denoising on single-frame, low-PSNR microscopy images even when the contrast is far out of distribution. While both low-contrast and high-contrast images show decreased denoising performance as input PSNR increases (Figs. 3 and 9), denoising improvement decreases at a faster rate for low-contrast images. This defines the out-of-distribution limit for the trained Noise2Noise plugin CNN: exceedingly low-contrast yet high-PSNR input images may have artifacts introduced into the image.

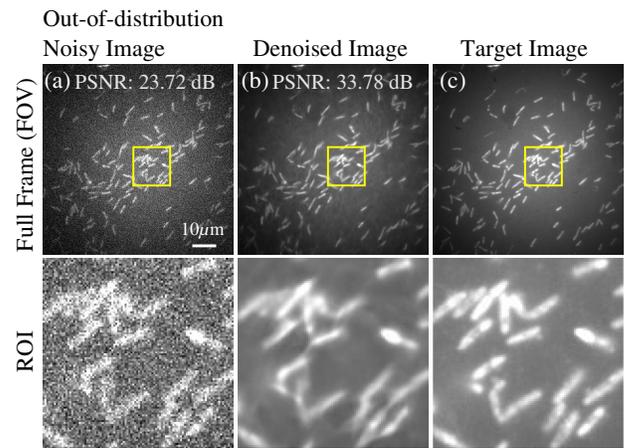

**Fig. 10.** Noise2Noise plugin denoising on *M. xanthus* autofluorescence images captured using a widefield fluorescence microscope. (a) Input noisy image, (b) image denoised using the Noise2Noise plugin, and (c) target image are represented by columns. The full-frame FOV is  $512 \times 512$  pixels; the selected ROI (yellow box) is  $150 \times 150$  pixels. Excitation wavelength 458 nm; emission wavelength 500–580 nm; sample excitation power 0.1 mW; integration time 100 ms; pixel size 150 nm.

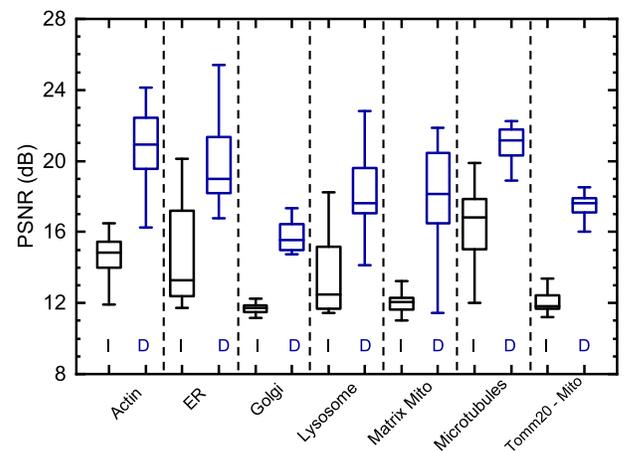

**Fig. 11.** PSNR distribution of noisy and denoised images taken from the out-of-distribution samples from the 3D RCAN dataset [33]. Gray and blue statistics indicate the box plot results of the input and denoised PSNR, respectively. I, input; D, denoised.

In summary, the Noise2Noise plugin demonstrates high-performance denoising on out-of-distribution low-PSNR and low-contrast ratio images taken from datasets external to the training dataset. By using bias-free convolutional layers, denoising on low-PSNR images is enabled. The trained CNN performs well on low-PSNR and low-contrast images outside of the dataset; however, the CNN introduces artifacts that degrade images in the *M. xanthus* dataset when contrast is low ( $<0.6$ ) and PSNR  $\gtrsim 36$  dB. However, such high-PSNR images would likely not be candidates for denoising, and therefore applying the Noise2Noise plugin is reasonable for practical experimental fluorescence microscopy PSNR and contrast. The Noise2Noise plugin shows exceptional performance ( $>8$  dB) for single-frame, low-PSNR images ( $<25$  dB) with typical experimental contrast (at least  $>0.5$ ).

### C. Out-of-Distribution Structures

While the ML model was trained on the FMD dataset (containing fluorescent-labeled BPAE cells and intravital genetically encoded fluorescent mouse brain and zebrafish samples), the previous section demonstrated that the model can perform accurate denoising on out-of-distribution *M. xanthus* structures (cellular features that do not appear in the training set images). To further evaluate performance on out-of-distribution structure, the denoising plugin performance was evaluated on the 3D RCAN dataset [33], consisting of fixed U2OS cells (cultivated from the human osteosarcoma bone tissue) transfected with mEmerald-Tomm20 labeling the outer mitochondrial membrane. The organelles are fluorescence stained with Alexa Fluor 488 Phalloidin (actin); ERmoxGFP (endoplasmic reticulum, ER); GalT-GFP (golgi); LAMP1-mApple (lysosomes); and Mouse- $\alpha$ -Tubulin primary antibody, Donkey- $\alpha$ -Mouse Biotin secondary antibody, and Alexa Fluor 488 Streptavidin (microtubules). In addition to a variety of out-of-distribution structures, the images were acquired using an out-of-distribution imaging modality: “instant structured illumination microscopy” (iSIM) [48]. These structures are chosen due to the diversity of features across length scales and are used for validation in other ML models [5]. Hence, the 3D RCAN dataset samples were chosen to represent the out-of-distribution structures. To evaluate plugin performance with varying SNRs in the RCAN dataset, synthetic Gaussian noise was added to the high-SNR images in the dataset.

Next, the Noise2Noise plugin is applied to the noisy images, and the quantitative improvements in PSNR and SSIM are computed. To account for differing image intensities (dynamic range of pixels values) from the acquired low-SNR and high-SNR samples, normalization is required before calculating the quantitative results. Images are percentile normalized using the method described in [17], and normalized images are denoted as  $\hat{X} = \text{Norm}(X, \text{percentile-low}, \text{percentile-high})$ , where  $X$  is the noisy image, second argument “percentile-low” indicates the lowest value of  $X$  [calculated as  $\text{percentile}(X, \text{percentile-low})$ ] and third argument “percentile-high” indicates the highest value of  $X$  [calculated as  $\text{percentile}(X, \text{percentile-high})$ ]. For noisy low-SNR/denoised images, percentile-low and percentile-high are set to one and 99, respectively. Similarly, for target (high-SNR) images, percentile-low and percentile-high are set to 0.1 and 99.9, respectively.

Table 3 provides a summary of average PSNR and SSIM improvements after denoising on the 3D RCAN dataset [33]. Average PSNR improvements for each class of structures range from 3.83 dB to 6.39 dB, and each of the 207 images of the dataset yields a positive increase in PSNR ( $\Delta\text{PSNR} > 0$ ) and SSIM ( $\Delta\text{SSIM} > 0$ ).

Figure 11 shows the Noise2Noise plugin performance in terms of PSNR distribution of noisy (gray color) and denoised (blue color) samples from the 3D RCAN dataset. From Table 3, our Noise2Noise plugin shows better image denoising performance on the out-of-distribution structures commonly arising in biomedical imaging. Qualitative results on the 3D RAN dataset, additional out-of-distribution datasets (such as fluorescence microscopy samples of BPAE cells, membrane samples [49]), and a fixed-pattern noise dataset captured using out-of-distribution modality total internal reflection fluorescence (TIRF) and lattice light-sheet microscopy (LLSM) microscope systems [5] are shown

**Table 3. Average PSNR Improvement ( $\Delta\text{PSNR} = \text{denoised image PSNR, input image PSNR}$ ) and Average SSIM Improvement ( $\Delta\text{SSIM} = \text{denoised image SSIM, input image SSIM}$ ) of the Noise2Noise Plugin When Applied to Out-of-distribution Samples of Fluorescence Microscopy Images from the 3D RCAN Dataset<sup>a</sup>**

Sample	Num. Samples	$\Delta\text{PSNR}$ (dB)	$\Delta\text{SSIM}$
Actin	21	6.39 (2.10)	0.59 (0.06)
ER	24	4.91 (1.35)	0.32 (0.09)
Golgi	34	3.83 (0.84)	0.21 (0.05)
Lysosome	35	4.99 (1.84)	0.33 (0.09)
Matrix Protein Mitochondria	36	6.10 (2.20)	0.33 (0.16)
Microtubule	32	4.38 (2.01)	0.39 (0.11)
Tomm20 Mitochondria	25	5.35 (0.79)	0.33 (0.06)

<sup>a</sup>Population standard deviations of improvements are given in parentheses.

in Supplement 1 Note S4 and in our GitHub repository under the Model validation folder [50].

### D. Application Specific Model Validation

The use of pre-trained neural network models bears the risk of introducing artifacts, especially when applied to images that differ from the training dataset. While this paper evaluates Noise2Noise plugin performance on some out-of-distribution contrast, SNR, structure, and imaging modality datasets, applying the plugin generally to untested data requires users to examine performance. The following steps can be applied to evaluate performance on new data sources, samples, or microscopy modalities. Application-specific validation follows three steps: image acquisition, applying denoising through the ML plugin, and performance verification.

First, obtain several single-channel low-SNR frames of the same fluorescence microscopy image FOV (i.e., at least greater than 10 2D sections at low excitation power or fast data acquisition rates). In addition to each individual frame, an estimate of the clean image is to be obtained either by acquiring the same FOV under long-term exposure or averaging all the acquired low-SNR frames together. To obtain the average, the collection of individual frames can be combined in Fiji/ImageJ as a “Stack,” then averaged using the Stack menu commands. Alternatively, a high SNR estimate can be obtained by performing averaging on the individual frames using a variety of software tools (e.g., Python, MATLAB). The output representing the clean image is then saved to a disk.

Second, load individual frames into Fiji/ImageJ and apply the Noise2Noise plugin. The installation of the Noise2Noise plugin and other details are provided in Supplement 1 Note S5 and the GitHub repository [50]. Save the output denoised image files to a disk.

Third, evaluate denoising performance by comparing the output of the denoising plugin to the estimated clean image. To evaluate performance qualitatively, first load both images simultaneously into Fiji/ImageJ while ensuring that both images are displayed with the same “Brightness/Contrast” scale (e.g., in the menu Image  $\rightarrow$  Adjust  $\rightarrow$  Brightness/Contrast...). Next, compare the structure of specific features of interest between the two images and examine the background (i.e., low-intensity) regions of the image for apparent artifacts. Quantitative verification can be performed by computing and comparing the PSNR and SSIM metrics of denoised images with respect to clean images. PSNR and SSIM for the raw images compared to the clean image and

the plugin output compared to the clean image can be calculated using ImageJ plugins [51,52]. The PSNR and SSIM of the output of the denoising plugin are expected to achieve at least  $\approx 3-8$  dB improvement in PSNR and a positive increase in SSIM compared to raw noisy images. Additional examples and detailed explanations of Noise2Noise plugin are provided in [50].

### E. Out-of-Distribution Data Dimensionality

Modern microscopy techniques commonly acquire data in four or more dimensions (three spatial, one temporal, and possibly more channels) and represent the data as 2D sections per channel at discrete acquisition times. It is therefore important to consider the model in the context of data with dimensionality outside of distribution from the training dataset. To address this, the model was trained on 2D sections of 3D spatial data of one to three channels, captured at specific times during dynamic intravital processes (e.g., blood flow and respiration). The noise of each individual channel can be described as MPG noise, and therefore, it is appropriate to apply the method to 2D sections of higher-dimension data using the MPG denoising method. Since the trained Noise2Noise plugin model performs denoising solely using information from a single 2D plane of the 3D sample, potential further performance improvements could be achieved by developing and training models using complete high-dimension data (e.g., the entire 3D stack in time [53]). Such an approach relies on significant computational improvements and would effectively amplify signals correlated in 3D space and time, leading to even improved denoising performance.

## 5. CONCLUSION

In this paper, we have presented an open-source software tool for fluorescence microscopy image denoising in real time (within tens of milliseconds) to denoise MPG noise. The new image denoising technique using the trained ML models (Noise2Noise plugin and DnCNN plugin) allows imaging eight times faster than the fundamental speed limit of fluorescence microscopy while maintaining the same PSNR. The system was trained using the experimentally obtained FMD dataset to develop a neural network that effectively removes MPG noise presented in widefield, confocal, and two-photon fluorescence microscopy images across a variety of cell types. We have shown here that the proposed Noise2Noise plugin ML model demonstrates larger PSNR improvements compared to the existing image denoising methods, performs denoising faster, and has been validated on a wide variety of images with varying SNRs, contrasts, structures, and imaging modalities.

**Funding.** Office of Science (0019312); Division of Chemical, Bioengineering, Environmental, and Transport Systems (1554516).

**Acknowledgment.** We thank Prof. Joshua D. ShROUT, Department of Civil and Environmental Engineering and Earth Sciences, University of Notre Dame, for providing *M. xanthus* cells for our experiments. Yide Zhang's research was supported by the Berry Family Foundation Graduate Fellowship of Advanced Diagnostics & Therapeutics (AD&T), University of Notre Dame. The authors further acknowledge Notre Dame's Center for Research Computing (CRC) for the use of GPUs to train the ML models and the Integrated Imaging Facility (NDIIF) for the use of the Nikon A1R MP confocal microscope and Nikon Eclipse 90i widefield microscope in NDIIF's Optical Microscopy Core.

**Disclosures.** The authors declare no conflicts of interest.

**Data availability.** Data and Code Availability. The FMD dataset mentioned in this paper is publicly available in [25]. The code for training the Noise2Noise plugin and DnCNN plugin architectures and underlying the results presented in the paper are publicly available in [43]. Also, this repository includes the estimation of noise parameters using MATLAB codes, Noise2Noise plugin validation on the W2S dataset, Noise2Noise plugin validation on out-of-distribution samples, and plugin source code in java with ImageJ plugins. Additional details of our Noise2Noise plugin denoising validation on the large datasets are available in [50]. The out-of-distribution structure dataset (3D RCAN dataset) used to generate underlying denoising results presented in this paper is available in [33].

**Supplemental document.** See [Supplement 1](#) for supporting content.

## REFERENCES

1. J. W. Lichtman and J.-A. Conchello, "Fluorescence microscopy," *Nat. Methods* **2**, 910–919 (2005).
2. E. H. Stelzer, "Contrast, resolution, pixelation, dynamic range and signal-to-noise ratio: fundamental limits to resolution in fluorescence light microscopy," *J. Microsc.* **189**, 15–24 (1998).
3. V. Mannam, Y. Zhang, X. Yuan, C. Ravasio, and S. S. Howard, "Machine learning for faster and smarter fluorescence lifetime imaging microscopy," *J. Phys. Photon.* **2**, 042005 (2020).
4. Y. Zhang, A. A. Khan, G. D. Vigil, and S. S. Howard, "Investigation of signal-to-noise ratio in frequency-domain multiphoton fluorescence lifetime imaging microscopy," *J. Opt. Soc. Am. A* **33**, B1–B11 (2016).
5. B. Mandracchia, X. Hua, C. Guo, J. Son, T. Urner, and S. Jia, "Fast and accurate sCMOS noise correction for fluorescence microscopy," *Nat. Commun.* **11**, 94 (2020).
6. K. Suhling, R. W. Airey, and B. L. Morgan, "Minimization of fixed pattern noise in photon event counting imaging," *Rev. Sci. Instrum.* **73**, 2917–2922 (2002).
7. T. J. Lambert and J. C. Waters, "Assessing camera performance for quantitative microscopy," *Methods Cell Biol.* **123**, 35–53 (2014).
8. A. Buades, B. Coll, and J.-M. Morel, "A non-local algorithm for image denoising," in *IEEE Computer Society Conference on Computer Vision and Pattern Recognition* (2005), Vol. 2, pp. 60–65.
9. K. Dabov, A. Foi, V. Katkovnik, and K. Egiazarian, "Image denoising by sparse 3-D transform-domain collaborative filtering," *IEEE Trans. Image Process.* **16**, 2080–2095 (2007).
10. M. Aharon, M. Elad, and A. Bruckstein, "K-SVD: an algorithm for designing overcomplete dictionaries for sparse representation," *IEEE Trans. Signal Process.* **54**, 4311–4322 (2006).
11. D. Zoran and Y. Weiss, "From learning models of natural image patches to whole image restoration," in *IEEE International Conference on Computer Vision* (2011), pp. 479–486.
12. S. Gu, L. Zhang, W. Zuo, and X. Feng, "Weighted nuclear norm minimization with application to image denoising," in *IEEE Conference on Computer Vision and Pattern Recognition* (2014), pp. 2862–2869.
13. M. Makitalo and A. Foi, "Optimal inversion of the generalized anisotropic transformation for Poisson-Gaussian noise," *IEEE Trans. Image Process.* **22**, 91–103 (2012).
14. C.-A. Deledalle, F. Tupin, and L. Denis, "Poisson NL means: unsupervised non local means for Poisson noise," in *IEEE International Conference on Image Processing* (2010), pp. 801–804.
15. J. Salmon, Z. Harmany, C.-A. Deledalle, and R. Willett, "Poisson noise reduction with non-local PCA," *J. Math. Imaging Vis.* **48**, 279–294 (2014).
16. F. Luisier, T. Blu, and M. Unser, "Image denoising in mixed Poisson-Gaussian noise," *IEEE Trans. Image Process.* **20**, 696–708 (2010).
17. M. Weigert, U. Schmidt, T. Boothe, *et al.*, "Content-aware image restoration: pushing the limits of fluorescence microscopy," *Nat. Methods* **15**, 1090–1097 (2018).
18. K. Zhang, W. Zuo, Y. Chen, D. Meng, and L. Zhang, "Beyond a Gaussian denoiser: residual learning of deep CNN for image denoising," *IEEE Trans. Image Process.* **26**, 3142–3155 (2017).
19. J. Lehtinen, J. Munkberg, J. Hasselgren, S. Laine, T. Karras, M. Aittala, and T. Aila, "Noise2Noise: learning image restoration without clean data," in *International Conference on Machine Learning* (PMLR, 2018), pp. 2965–2974.

20. A. Krull, T.-O. Buchholz, and F. Jug, "Noise2Void-learning denoising from single noisy images," in *IEEE/CVF Conference on Computer Vision and Pattern Recognition* (2019), pp. 2129–2137.
21. A. Krull, T. Vičar, M. Prakash, M. Lalit, and F. Jug, "Probabilistic Noise2Void: unsupervised content-aware denoising," *Front. Comput. Sci.* **2**, 5 (2020).
22. J. Lequyer, R. Philip, A. Sharma, and L. Pelletier, "Noise2Fast: fast self-supervised single image blind denoising," arXiv:2108.10209 (2021).
23. W. Khademi, S. Rao, C. Minnerath, G. Hagen, and J. Ventura, "Self-supervised Poisson-Gaussian denoising," in *IEEE/CVF Winter Conference on Applications of Computer Vision* (2021), pp. 2131–2139.
24. A. S. Goncharova, A. Honigmann, F. Jug, and A. Krull, "Improving blind spot denoising for microscopy," in *European Conference on Computer Vision* (Springer, 2020), pp. 380–393.
25. V. Mannam, Y. Zhang, Y. Zhu, and S. S. Howard, "Fluorescence microscopy denoising (FMD) dataset," CurateND (2019), <https://doi.org/doi:10.7274/r0-ed2r-4052>.
26. Y. Zhang, Y. Zhu, E. Nichols, Q. Wang, S. Zhang, C. Smith, and S. Howard, "A Poisson-Gaussian denoising dataset with real fluorescence microscopy images," in *IEEE Conference on Computer Vision and Pattern Recognition* (2019), pp. 11710–11718.
27. R. Zhou, M. El Helou, D. Sage, T. Laroche, A. Seitz, and S. Süsstrunk, "W2S: microscopy data with joint denoising and super-resolution for widefield to SIM mapping," in *European Conference on Computer Vision* (Springer, 2020), pp. 474–491.
28. J. Pawley, *Handbook of Biological Confocal Microscopy* (Springer, 2006), Vol. **236**.
29. W. Denk, J. H. Strickler, and W. W. Webb, "Two-photon laser scanning fluorescence microscopy," *Science* **248**, 73–76 (1990).
30. P. J. Verveer, M. J. Gemkow, and T. M. Jovin, "A comparison of image restoration approaches applied to three-dimensional confocal and wide-field fluorescence microscopy," *J. Microsc.* **193**, 50–61 (1999).
31. H. Wang, Y. Rivenson, Y. Jin, Z. Wei, R. Gao, H. Günaydn, L. A. Bentolila, C. Kural, and A. Ozcan, "Deep learning enables cross-modality super-resolution in fluorescence microscopy," *Nat. Methods* **16**, 103–110 (2019).
32. W. Ouyang, A. Aristov, M. Lelek, X. Hao, and C. Zimmer, "Deep learning massively accelerates super-resolution localization microscopy," *Nat. Biotechnol.* **36**, 460–468 (2018).
33. J. Chen, H. Sasaki, H. Lai, Y. Su, J. Liu, Y. Wu, A. Zhovmer, C. A. Combs, I. Rey-Suarez, H.-Y. Chang, C. C. Huang, X. Li, M. Guo, S. Nizambad, A. Upadhyaya, S.-J. J. Lee, L. A. G. Lucas, and H. Shroff, "Three-dimensional residual channel attention networks denoise and sharpen fluorescence microscopy image volumes," *Nat. Methods* **18**, 678–687 (2021).
34. S. Fujita and X.-H. Han, "Cell detection and segmentation in microscopy images with improved mask R-CNN," in *ACCV Workshops* (2020), pp. 58–70.
35. J. Schindelin, C. T. Rueden, M. C. Hiner, and K. W. Eliceiri, "The ImageJ ecosystem: An open platform for biomedical image analysis," *Mol. Reprod. Dev.* **82**, 518–529 (2015).
36. C. Rueden, D. Schmidt, and B. Wilhelm, "ImageJ Tensorflow library," 2017, <https://github.com/imagej/imagej-tensorflow>.
37. J. Boulanger, C. Kervrann, P. Bouthemy, P. Elbau, J.-B. Sibarita, and J. Salamero, "Patch-based nonlocal functional for denoising fluorescence microscopy image sequences," *IEEE Trans. Med. Imaging* **29**, 442–454 (2009).
38. E. Gómez-de Mariscal, C. Garca-López-de Haro, L. Donati, M. Unser, A. Muñoz-Barrutia, and D. Sage, "DeepImageJ: a user-friendly plugin to run deep learning models in ImageJ," *Nat. Methods* **18**, 1192–1195 (2021).
39. L. von Chamier, R. F. Laine, J. Jukkala, et al., "Democratising deep learning for microscopy with ZeroCostDL4Mic," *Nat. Commun.* **12**, 2276 (2021).
40. A. Foi, M. Trimeche, V. Katkovnik, and K. Egiazarian, "Practical Poissonian-Gaussian noise modeling and fitting for single-image raw-data," *IEEE Trans. Image Process.* **17**, 1737–1754 (2008).
41. Z. Wang, A. C. Bovik, H. R. Sheikh, and E. P. Simoncelli, "Image quality assessment: from error visibility to structural similarity," *IEEE Trans. Image Process.* **13**, 600–612 (2004).
42. A. Hore and D. Ziou, "Image quality metrics: PSNR vs. SSIM," in *20th International Conference on Pattern Recognition* (IEEE, 2010), pp. 2366–2369.
43. <https://github.com/ND-HowardGroup/Instant-Image-Denoising>.
44. A. Krull, T.-O. Buchholz, and F. Jug, "ImageJ Noise2Void image denoising," 2018, <https://imagej.net/plugins/n2v>.
45. S. Mohan, Z. Kadkhodaie, E. P. Simoncelli, and C. Fernandez-Granda, "Robust and interpretable blind image denoising via bias-free convolutional neural networks," in *International Conference on Learning Representations* (ICLR, 2019).
46. V. Sundaresan, A. R. Cutri, J. Metro, C. S. Madukoma, J. D. ShROUT, A. J. Hoffman, K. A. Willets, and P. W. Bohn, "Potential dependent spectroelectrochemistry of electrofluorogenic dyes on indium-tin oxide," *Electrochem. Sci. Adv.* **2**, e2100094 (2021).
47. C. W. Harvey, C. S. Madukoma, S. Mahserejian, M. S. Alber, and J. D. ShROUT, "Cell division resets polarity and motility for the bacterium *Myxococcus xanthus*," *J. Bacteriol.* **196**, 3853–3861 (2014).
48. A. G. York, P. Chandris, D. Dalle Nogare, J. Head, P. Wawrzusin, R. S. Fischer, A. Chitnis, and H. Shroff, "Instant super-resolution imaging in live cells and embryos via analog image processing," *Nat. Methods* **10**, 1122–1126 (2013).
49. G. M. Hagen, J. Bendesky, R. Machado, T.-A. Nguyen, T. Kumar, and J. Ventura, "Fluorescence microscopy datasets for training deep neural networks," *GigaScience* **10**, giab032 (2021).
50. [https://github.com/ND-HowardGroup/Instant-Image-Denoising/tree/master/Plugins/Model\\_validation](https://github.com/ND-HowardGroup/Instant-Image-Denoising/tree/master/Plugins/Model_validation).
51. D. Sage, "ImageJ's plugin to assess the quality of images," 2017, <http://bigwww.epfl.ch/sage/soft/snr/>.
52. P. R. Gabriel, "SSIM: a Java plugin in ImageJ," 2008, <https://imagej.nih.gov/ij/plugins/ssim-index.html>.
53. M. Papkov, K. Roberts, L. A. Madisson, J. Shilts, O. Bayraktar, D. Fishman, K. Palo, and L. Parts, "Noise2Stack: improving image restoration by learning from volumetric data," in *International Workshop on Machine Learning for Medical Image Reconstruction* (Springer, 2021), pp. 99–108.