

Video Article

Live Confocal Imaging of Developing Arabidopsis Flowers

Nathanaël Prunet¹

¹Department of Biology and Biological Engineering, California Institute of Technology

Correspondence to: Nathanaël Prunet at nprunet@caltech.edu

URL: <https://www.jove.com/video/55156>

DOI: [doi:10.3791/55156](https://doi.org/10.3791/55156)

Keywords: Developmental Biology, Issue 122, live confocal imaging, confocal microscopy, Arabidopsis, flower, flower development, flower meristem

Date Published: 4/1/2017

Citation: Prunet, N. Live Confocal Imaging of Developing Arabidopsis Flowers. *J. Vis. Exp.* (122), e55156, doi:10.3791/55156 (2017).

Abstract

The study of plant growth and development has long relied on experimental techniques using dead, fixed tissues and lacking proper cellular resolution. Recent advances in confocal microscopy, combined with the development of numerous fluorophores, have overcome these issues and opened the possibility to study the expression of several genes simultaneously, with a good cellular resolution, in live samples. Live confocal imaging provides plant biologists with a powerful tool to study development, and has been extensively used to study root growth and the formation of lateral organs on the flanks of the shoot apical meristem. However, it has not been widely applied to the study of flower development, in part due to challenges that are specific to imaging flowers, such as the sepals that grow over the flower meristem, and filter out the fluorescence from underlying tissues. Here, we present a detailed protocol to perform live confocal imaging on live, developing Arabidopsis flower buds, using either an upright or an inverted microscope.

Video Link

The video component of this article can be found at <https://www.jove.com/video/55156/>

Introduction

Most of the plant body forms post-embryonically from groups of stem cells situated at or near the tip — or apical meristems — of the shoots and roots. All the above-ground structures of the adult plant derive from the shoot apical meristem (SAM), which continuously produces lateral organs on its flanks: leaves during the vegetative phase, and flower meristems (FMs) after it transitions to the reproductive phase. FMs in turn develop into flowers. While the Arabidopsis SAM produces lateral organs one at a time, in an iterative, spiral pattern, FMs produce four types of floral organs within four whorls, in a partially synchronous manner, with multiple developmental programs unraveling simultaneously. The genetic networks underlying the specification of the identity of the different floral organs have been partially deciphered (for reviews, see references^{1,2}), but many aspects of flower development, such as floral organ positioning and the definition of boundaries between whorls, remain poorly understood.

Early molecular genetic studies of plant development mostly relied on techniques such as *in situ* hybridization and GUS reporters to analyze gene expression. While these methods have provided a wealth of information and greatly contributed to our understanding of plant growth and flower development, there are important limitations: they lack good cellular resolution, do not allow for the easy observation of the expression pattern of multiple genes in the same samples, and importantly, can only be applied to dead, fixed tissues. Recent advances in confocal microscopy have overcome these limitations, and provide developmental biologists with a powerful tool to investigate the processes underlying plant morphogenesis. In particular, confocal microscopy allows for the observation of live tissues and organs throughout their formation, which is critical to fully understand a quintessentially dynamic process such as development.

Live confocal imaging has been extensively used to analyze the aerial growth of plants and the production of lateral organs by the SAM (e.g. references^{3,4,5,6}), but with the exception of a few reports (e.g. references^{7,8,9,10}), it has not been widely applied to the study of flower development. Protocols for live confocal imaging of the SAM are available (e.g. references^{11,12}), and provide a good base for how to image the developing flower buds that surround the SAM. However, imaging flower buds presents specific challenges: for instance, flower buds quickly become bigger than the SAM, and at stage 4, sepals start covering the FM (stages as described in reference¹³), and dimming the fluorescence from underlying tissues (**Figure 2A**). Here, we provide a detailed protocol explaining how to perform live confocal imaging on developing Arabidopsis flower buds using either an upright or an inverted microscope and a water-dipping lens. We previously published this protocol in reference¹⁴.

Protocol

1. Media and Dishes Preparation

1. Prepare dissecting dishes by filling round plastic boxes (approximately 6 cm wide, 2 cm deep) to 0.5 cm with 2% agarose.
2. **Prepare imaging dishes.**

1. For an upright confocal microscope, fill rectangular plastic hinged boxes (approximately 7 cm long, 4.5 cm wide, 3 cm deep) to 0.5 cm with imaging medium (see section 1.3, "Imaging medium"; **Figure 1F**).
2. For an inverted confocal microscope: fill small Petri dish (approximately 3.5 cm wide, 1 cm deep) exactly to the brim with imaging medium (see section 1.3, "Imaging medium"; **Figure 1G**).
3. **Prepare imaging medium.**
 1. For single point imaging, use 1% agarose.
 2. For time-lapse experiments, use apex growth medium (0.5x Murashige and Skoog basal salt mixture without vitamin, 1% sucrose, 0.8% agarose, pH 5.8 with potassium hydroxide solution, supplemented with vitamin [0.01% myo-inositol, 0.0001% nicotinic acid, 0.0001% pyridoxine hydrochloride, 0.001% thiamine hydrochloride, 0.0002% glycine] and cytokinin [500 nM N6-benzyladenine])¹⁵.

2. Plant Growth

1. Sow sterilized seeds on MS plates (0.5 or 1x Murashige and Skoog basal salt mixture without vitamin, 0.8% agar, pH 5.8 with potassium hydroxide solution) with appropriate selection. Place plates in long day (16 h light) or continuous day, 16-22 °C conditions for two weeks.
2. Transplant seedlings onto soil with sufficient spacing to allow robust development. Place plants in short day, 16-22 °C conditions for three weeks.
NOTE: Growing plants in long day or continuous day conditions directly after transplanting results in premature flowering and inflorescences that are less vigorous, and much harder to dissect. Similarly, leaving plants under short day conditions for more than three weeks results in day length-independent flowering and inflorescences that are less vigorous.
3. Transfer plants to long day (16 h light) or continuous day, 16-22 °C conditions until they flower. Shoot apices are easiest to dissect when the inflorescence is 2-10 cm long.

3. Dissection of the Shoot Apex

1. Optional: Using a fine sharpening stone and a drop of light oil, sharpen forceps under the stereomicroscope to make them blade-like, not point-like.
2. Remove siliques, older flower buds and secondary inflorescences from the primary inflorescence by pushing the base of the peduncles with the forceps until they break. Remove as many flowers as possible without magnification (**Figure 1**, compare A and B).
3. Using forceps, pierce a vertical hole in the agarose of a dissecting dish, cut off the last 0.5 cm of the inflorescence and stick it vertically in the agarose. The remaining flower buds must be above the agarose surface.
4. Fill the imaging dish with sterile, de-ionized water so that the shoot apex is fully immersed. Place the dissecting dish under the stereomicroscope, and remove the air trapped around the shoot apex by creating water jets with a 1,000 µL pipette (**Figure 1**, compare C and D).
5. Under the stereomicroscope, use the forceps to remove flower buds that are not to be imaged by pushing on the base of the peduncle or the top of the bud until the peduncle breaks (**Figure 1E**). It is important that the peduncles break cleanly at the junction with the stem, as leftover peduncles hinder the removal of the young flower buds.
6. If imaging only flower buds stage 5 and younger (**Figure 2A**), remove the water before dissecting stage 6-8 flower buds. If imaging flower buds stage 5 and older, proceed to sepal ablation (see section 5.1). Otherwise, proceed to step 3.7.
7. Using forceps, pierce a vertical hole in the medium of an imaging dish, and stick the dissected apex upright in the medium, so that only the shoot apical meristem and surrounding flower buds are above the surface of the medium. Depending on which flower bud(s) are to be imaged, it may be necessary to slightly tilt the sample, as flower buds are not necessarily oriented exactly like the stem (**Figure 2A**).
8. Proceed to staining the sample if needed. If not staining and/or imaging the sample right away, add water and close the imaging dish to prevent dehydration. If using an inverted microscope, place the imaging dish in a transparent plastic box and close it.

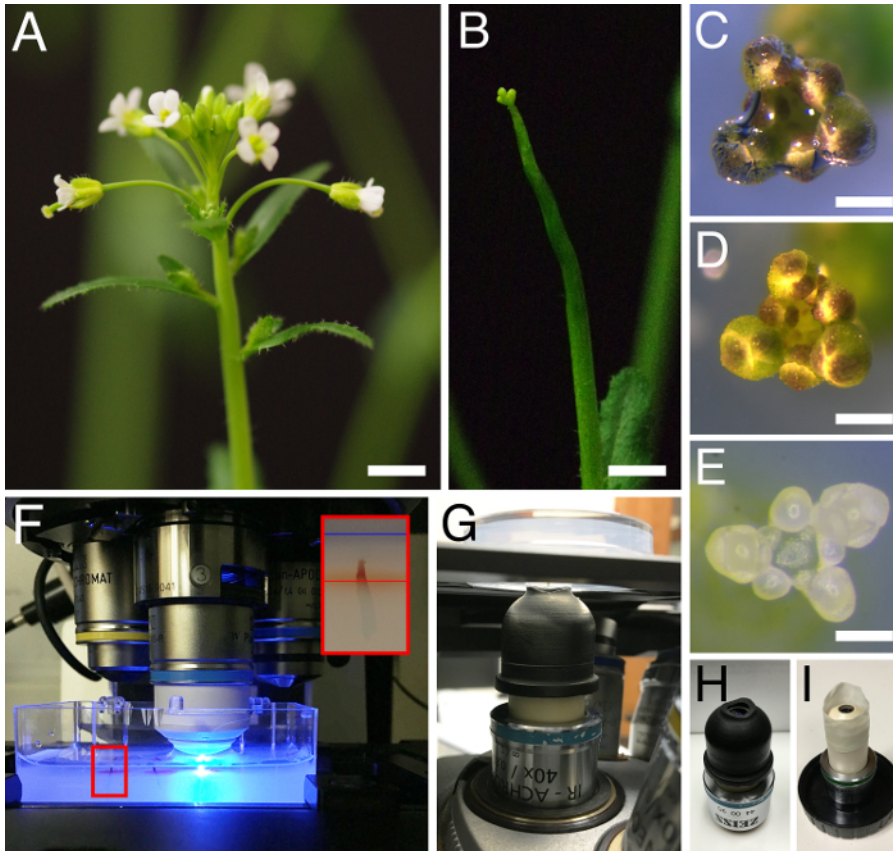


Figure 1: Preparation of the shoot apices for live confocal imaging. (A-E) Dissection of a shoot apex for imaging. Inflorescence before (A) and after (B) removal of silicles and older flowers. (C-D) Shoot apex immersed in water in the dissecting dish, with an air bubble trapped at the tip (C) and after removal of the air bubble (D). (E) Shoot apex in the dissecting dish, after dissection of flower buds older than stage 5. (F-G) View of shoot apices in the imaging dish, on the stage of an upright (F) and inverted (G) confocal microscope. In (F), a 40X water-dipping lens is positioned above one of the apices, with the tip of the lens immersed in water. In (G), the shoot apex is positioned upside down above the 40X water-dipping lens, with a water column connecting the imaging medium to the tip of the lens. The smaller panel in (F) shows a higher magnification view of the area in the red rectangle, with a shoot apex inserted in imaging medium; red and blue lines indicate the surface of the medium and water, respectively. (H-I) Examples of custom-made devices that allow adding more water at the tip of the water-dipping lens: a silicon rubber sleeve made from a spark plug boot (H), and a makeshift sleeve made from the finger of a powderless latex glove (G). Scale bars = 0.5 cm in (A) and (B), 0.1 cm in (C) and (D), and 100 μ m in (E). This figure was initially published in reference¹⁴. [Please click here to view a larger version of this figure.](#)

4. Staining

NOTE: Cell walls or plasma membranes can be stained with propidium iodide or FM4-64, respectively, to achieve cellular resolution during imaging. Alternatively, reporter lines with fluorescent proteins tagged to the plasma membrane can be used^{6,16}.

1. Remove water from the imaging dish, and ensure the surface of the medium is dry to avoid diluting the dye.
2. Under the stereomicroscope, apply 20 to 30 μ L of either 1 mg/mL propidium iodide solution, or 80 μ g/mL FM4-64 solution to the dissected apex with a 10 μ L pipette. Make sure the whole sample is covered in dye to ensure a homogeneous staining.
3. Stain for 2 min if using propidium iodide, or 20 min if using FM4-64.
4. Rinse twice with sterile, deionized water.

5. Sepal Ablation

NOTE: At stage 4, sepal primordia start covering the FM. As they grow, they filter out the fluorescence from underlying tissue, and hinder the imaging process. Sepals can either be removed manually using a metal pin mounted on a pin-vise, or prevented to grow using laser ablation on emerging sepal primordia. Alternatively, it is possible in some cases to use mutants such as *apetala1-1*, in which sepals are missing or replaced by leaf-like organs that do not cover the FM while the central part of the flower develops normally (Figure 2F)¹⁷.

1. **Sepal ablation on flower buds stage 5 and older using a pin-vise holding a straight metal needle (Figure 2, E1-E2)**
 1. Place the dissecting dish with the dissected apex under the stereomicroscope, and set the magnification to maximum. Immerse the shoot apex in sterile, deionized water, or alternatively, use a 1,000 μ L pipette to apply water to the shoot apex regularly while dissecting the sepals.

2. With the pin-vise, position the pin on top of the abaxial sepal, tangentially relative to the shoot apex. Gently push the sepal away from the shoot apex until it breaks away from the flower bud.
3. Proceed similarly with the adaxial sepal, but push the sepal towards the shoot apex.
4. Proceed similarly with the lateral sepals, but position the pin radially relative to the shoot apex, and push the sepals to the side, away from the flower bud.
5. Immerse the shoot apex in sterile, deionized water for a few minutes to prevent dehydration.

2. Laser ablation of sepal primordia at stage 3-4 (Figure 2, C5-D5)

1. Place the imaging dish with the stained, dissected apex on the confocal microscope stage. Locate and focus on the apex as if preparing to image it (see section 6, "imaging setup").
2. Using the laser ablation system software, define the ablation zone, which corresponds to the crest and tip of the emerging sepal primordia before they start covering the flower meristem (typically, at stage 3 for abaxial and adaxial sepals, and at late stage 3/early stage 4 for lateral sepals).
3. Set laser power and dwelling time to appropriate parameters to ablate enough cells without inflicting too much damage to the underlying tissues. Typically, depending on the laser ablation system used, proper parameters initially need to be identified through a trial-and-error process to ensure that enough cells are removed to prevent the sepals to subsequently grow over the FM, without affecting the rest of the flower bud. Using too much laser power and dwelling time results in damages to the center of the flower, affecting its growth and/or survival. Once the proper parameters have been identified, they can be reused for subsequent experiments. NOTE: If the initial ablation proves insufficient, it is possible to proceed to a second ablation on the following days (Figure 2, C4 and D2).

6. Imaging Setup

1. Use an upright microscope.

1. Fill the imaging dish with the dissected apices with sterile, deionized water so that the surface of the medium is covered in 2.5-5 mm water (Figure 1F). Completely immerse the samples.
NOTE: When imaging with an upright microscope, several apices can be placed in the same imaging dish. However, if the imaging process lasts more than an hour, propidium iodide tends to be diluted, and some samples might need re-staining prior to imaging.
2. Place the imaging dish on the microscope stage (Figure 1F). Lower the water-dipping lens and raise the stage so that the tip of the lens dips in the water. Proceed cautiously so as not to crush the dissected apices or dip the lens into the imaging medium.
3. If an air bubble is trapped at the tip of the lens, make water jets with a 1,000 μ L pipette to remove it.
4. Under epifluorescence illumination, position one of the dissected apices into the lens field using the XY controller. Look through the eyepieces and focus on the sample using the Z controller. Proceed cautiously in order not to crush the sample.
5. Using the confocal microscope software, zoom on the flower bud to be imaged, and proceed to imaging your sample.

2. Use an inverted microscope.

1. Using a 1,000 μ L pipette, put a drop of sterile, deionized water on the tip of the lens.
2. Hold the imaging dish upside-down, and with a 1,000 μ L pipette, add a drop of sterile, deionized water to the dissected apex.
3. Place the imaging dish upside-down on the microscope stage (Figure 1G). Proceed cautiously so as not to crush the sample.
4. Under epifluorescence illumination, position the dissected apex over the tip of the lens using the XY controller. With the Z controller, carefully lower the stage until the sample reaches the drop of water at the tip of the lens. A water column should form between the tip of the lens and the medium.
5. If the water column does not form, carefully add a drop of water to the sample with a 1,000 μ L pipette. Add makeshift sleeves to the lens in order to add more water at its tip, which facilitates the establishment and maintenance of the water column (Figure 1H and 1I).
6. Under epifluorescence illumination, look through the eyepieces and focus on the sample using the Z controller. Proceed cautiously so as not to crush the sample.
7. Using the confocal microscope software, zoom on the flower bud to be imaged, and proceed to imaging the sample.

3. If performing time-lapse experiments, pour the water out of the imaging dish and close it to prevent dehydration in between time points. Place the imaging dish with the samples in long day (16 h light) or continuous day, 16-22 °C conditions. Re-stain samples before each time point. NOTE: While cells in developing floral organs can divide faster, cells in the floral meristem divide once or twice every second day, and cell lineages are easy to track with time points every 24 h. However, if need be, it is possible to image the samples every six hours.

7. Considerations on the Imaging Parameters

NOTE: How to set up the imaging parameters depends a lot on the confocal system used. Below are suggestions for some of these parameters that can be used with any confocal microscope. For more considerations on the imaging parameters, see the **Discussion** section and reference¹⁴.

1. For the XY resolution: use 1,024 x 1,024 for a good cellular resolution. Alternatively, use 512 x 512 to reduce imaging time, at the expense of the resolution.
2. For the Z resolution: set the step-size to 0.5-1.5 μ m. Lowering the step-size increases the Z resolution but also the imaging time.
3. For the pinhole: set the pinhole to 1-1.5 airy units. Increasing the pinhole increases the intensity of the signal, but also the noise from non-focal planes.

8. Visualization of Confocal Data

1. Turn still images of time points of flower development into a movie.

1. Open the still images (in jpeg or tif format) in the software (e.g., Fiji).
2. Go to Image > Stacks > Images to Stack. The open images will be grouped in a stack.
3. If the order of the images in the stack does not correspond to the chronological order, use the Stack Sorter plugin (Plugins > Stacks > Stack Sorter) to reorder them.
4. To turn the stack into a movie, go to File > Save As > AVI.

Representative Results

Figure 2 presents different views of confocal Z-stacks of live *Arabidopsis* flower buds expressing different fluorescent reporter genes, and stained with either propidium iodide (**Figure 2, A-C5 and G**) or FM4-64 (**Figure 2, E1-F**) to provide a clear cellular resolution. Most confocal systems allow for the imaging of two fluorophores with non-overlapping emission spectra such as GFP or YFP together with either propidium iodide or FM4-64 (**Figure 2A-2F**). The best confocal systems are also able to separate multiple fluorophores with partially overlapping emission spectra in the same samples, such as GFP and YFP, and dsRed and propidium iodide (**Figure 2G**).

Three examples of time-lapse experiments are presented in **Figure 2** and show flower buds developing normally after sepal ablation was performed either with a laser ablation system (**Figure 2, C1-C5 and D1-D5**) or manually with a pin-vise and a metal pin (**Figure 2, E1-E2**). Note that the flower bud shown in **Figure 2, E1-E2** is from a *superman-1* mutant plant, which is why carpel primordia are not observed at the center of the bud at stage 7. Laser ablation of emerging sepal primordia of the flower bud shown in **Figure 2, C1-C5** prevents them from covering the FM, which can still be imaged at stage 5 (**Figure 2, C5**). Conversely, in the case of the flower bud shown in **Figure 2, D1-D5**, laser ablation only delayed sepal growth, but sepals eventually grew to cover most of the FM by stage 5 (**Figure 2, D5**). On the contrary, if too much laser power and dwelling time is used during laser ablation, damage spreads to the central part of the flower and affects its growth and sometimes results in the subsequent death of the whole flower bud.

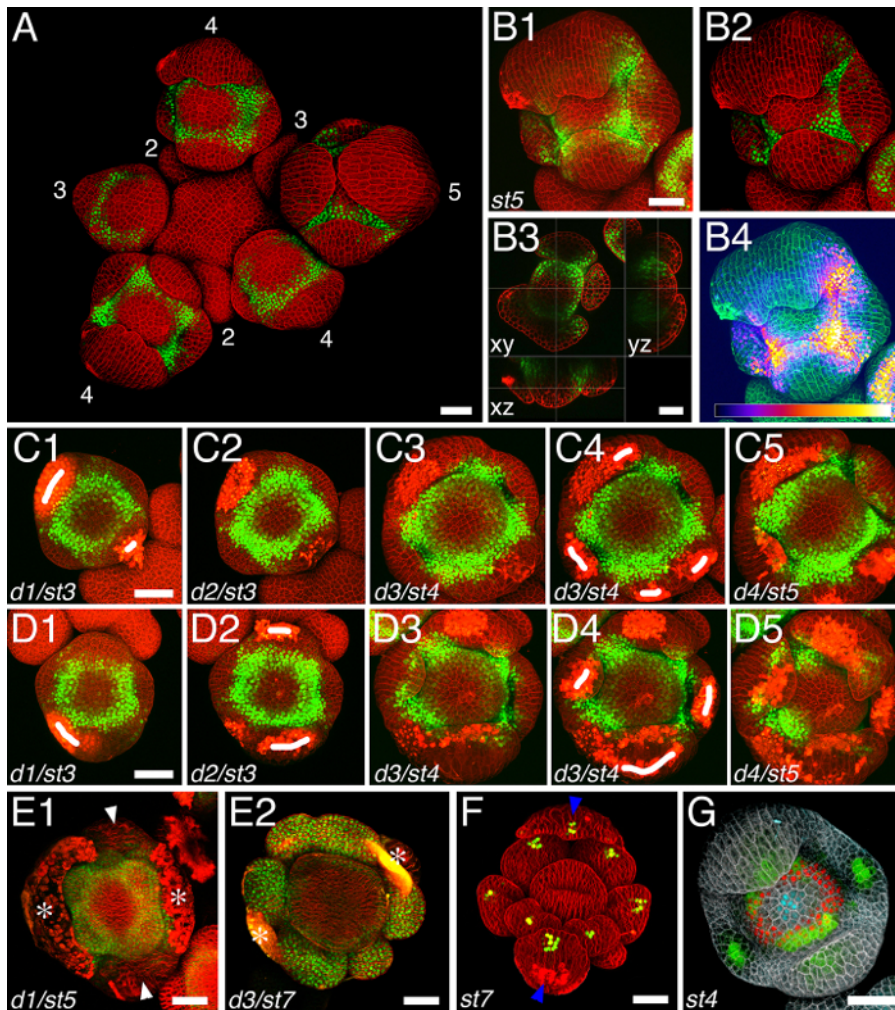


Figure 2: Visualization of confocal Z-stacks of live flower buds. A and B2 are transparency views; B1 and B4-G are maximum intensity projections; B3 is an orthogonal slice view. (A-E2) flowers expressing a Venus reporter for the *APETALA3* gene (green); cell walls were stained with propidium iodide (red, except for B4, green). (A) Inflorescence; numbers indicate floral stages; sepals in stage 4 and 5 flowers filter out the fluorescence of the Venus reporter, which normally forms a ring. Note that some flower buds appear tilted compared to the inflorescence. (B1-B4) Four different views of the same confocal Z-stack of a stage 5 flower bud: maximum intensity projections (B1 and B4) with Venus signal intensity coded with pixel intensity in the green color (B1) or with fire color code (B4; color code is shown at the bottom of the panel); transparency view (B2) and orthogonal slice view (B3; the XY, XZ and YZ orientation of the slices is indicated in the panel). (C1-C5) 4-day time-lapse of an individual flower bud from stage 3 to stage 5; laser ablations (marked as white traits) performed on day 1 and day 3 were sufficient to prevent the sepals from covering the center of the flower bud at stage 5. (D1-D5) 4-day time-lapse of an individual flower bud from stage 3 to stage 5; laser ablations performed on day 1, 2 and 3 were insufficient to prevent the sepals from covering the center of the flower bud. (E1-E2) Individual flower bud after manual removal of the abaxial and adaxial sepals (E1), and of all sepals (E2); white arrowheads indicate remaining sepals; white asterisks indicate scars resulting from the removal of the sepals. (F) stage 7 *apetala1-1* flower expressing a *DR5-3xVenusN7* reporter¹⁸ (green); plasma membranes were stained with FM4-64 (red); blue arrowheads indicate the leaf-like structures that replace sepals and do not cover the flower bud. (G) Stage 4 flower bud expressing a fluorescent GFP reporter for *DORNROSCHEN-LIKE1*⁷ (green), a Venus reporter for *SUPERMAN* (red) and a dsRed reporter for *CLAVATA3*¹⁹ (cyan); cells walls were stained with propidium iodide (grey). d: day; st: stage. Scale bars = 25 μ m; scale is identical in B1, B2 and B4, in C1-C5 and in D1-D5. This figure was modified from reference¹⁴. [Please click here to view a larger version of this figure.](#)

Discussion

Here, we provide a protocol for the imaging of live, developing *Arabidopsis* flower buds with a confocal microscope and a water-dipping lens. While this can be done with either an upright or an inverted microscope, it is easier and faster to use the former. It is also worth noting that different confocal systems vary significantly in terms of speed, sensitivity, and ability to separate wavelengths. The best confocal microscopes now allow to image several different channels (e.g. GFP, YFP, dsRed and propidium iodide; **Figure 2G**) in the same samples. Some confocal microscopes are equipped with a spectral detector, which can collect the fluorescence from several fluorophores simultaneously and separate them afterwards. When using a regular detector, however, a set of lasers and filters must be used to excite and separate the fluorescence from different fluorophores. Some fluorophores have fully distinct emission spectra (e.g. CFP and YFP), and can be imaged simultaneously. Other fluorophores have partially overlapping emission spectra (e.g. GFP and YFP), and usually need to be imaged separately, which increases imaging time. Imaging close fluorophores also often requires restricting how much of the spectrum is collected for each channel to prevent

fluorescence from one fluorophore from leaking into another channel. This causes a loss of intensity of the collected signal, which can be compensated by an increase in laser power and gain. However, prolonged imaging time and increased laser power may bleach and/or damage the sample. Increased laser power and/or gain may also increase background noise. Optimization of signal intensity, resolution and imaging time for each sample is done through a trial-and-error process, and involves permanent trade-offs.

This protocol explains how to perform live confocal imaging of flower buds developing on the flanks of a dissected shoot apex. It may be argued that the fact that the shoot apex is not connected to the rest of the plant and grows in a medium containing cytokinins might affect SAM and flower development. However, this medium was designed empirically through a trial-and-error process to ensure the dissected shoot apical meristems produce new flower buds at the normal plastochron, and that these flower buds develop normally¹⁵. Similarly, sepal dissection does not appear to affect gene expression patterns or the development of the rest of the flower. Other protocols detail how to perform live confocal imaging of the SAM still attached to the rest of the plant (e.g. reference¹¹). Such protocols could be adapted to, and offer a useful alternative for the imaging of developing flowers, particularly when studying cytokinin-related processes. They present several disadvantages however: the stem elongates and twists, and changes the orientation of the flower buds; and while imaging a shoot apex attached to the rest of the plant works well with an upright microscope, it is much more complicated to do so with an inverted microscope.

Immersion lenses have a better numerical aperture (NA) than "dry" lenses (i.e. lenses that are separated from the sample by air), and therefore provide a finer resolution of the sample, which is critical when using a confocal microscope. Using a water-dipping lens thus provides a much better NA than a dry lens, while preventing the shoot apex from dehydrating during the imaging process. While glycerin and oil lenses have an even higher NA than water lenses, they require the use of a coverslip, which is very impractical with a sample the size of a shoot apex, which also keeps growing during time-lapse experiments. Given the size of the samples (depending on the floral stages, stacks can be over 150 µm thick), it is also important to use a lens with a long working distance. We typically use a 20 or 40X lens with a NA of 1 and a working distance of 1.7-2.5 mm.

Confocal microscope software offers several ways to visualize confocal data, including three-dimensional reconstructions (**Figure 2, B1, B2 and B4**) and slice views (**Figure 2, B3**). While the most commonly used three-dimensional views are maximum intensity projections of the confocal Z-stacks (**Figure 2, B1 and B4**), some software (e.g. ZEN and Imaris) also offer transparency views that filters out the signal from the deeper parts of the sample (**Figure 2, B2**), which can be useful when looking at processes taking place, or genes expressed in, the epidermal layer. Signal intensity can either be coded with a single color (**Figure 2, B1**), using pixel intensity, or using a color code (**Figure 2, B4**).

Live confocal imaging not only offers valuable qualitative insights into flower development, it also potentially provides developmental biologists with a wealth of quantitative data. Different software (e.g. Imaris, Fiji, MARS-ALT, MorphoGraphX^{15,20,21}) allows for the quantification of the number or the volume of cells expressing one or several reporters, or the level of expression of a reporter in different cells. Such software can also be used to perform automatic cell segmentation, track cell lineages and quantify growth. This different software offers similar tools, but also differs in many ways. MARS-ALT and MorphoGraphX were designed specifically for plants^{15,20,21}, unlike Imaris and Fiji. MorphoGraphX only allows for the segmentation and analysis of the epidermal layer, while Imaris and MARS-ALT allow for the segmentation and analysis of the whole sample. However, MARS-ALT requires the prior imaging of the same sample from three different angles¹⁵, and Imaris only requires a single stack. Access to quantitative information is critical to further our understanding of flower development.

Disclosures

The author has nothing to disclose.

Acknowledgements

The author wish to thank Prof. Elliot M. Meyerowitz for his support, and comments on the manuscript, as well as Ann Lavanway at Dartmouth College and Dr. Andres Collazo at the Biological Imaging Facility at Caltech for their help in solving technical issues with the live confocal imaging. Nathanaël Prunet's work is supported by the US National Institutes of Health through Grant R01 GM104244 to Elliot M. Meyerowitz.

References

1. Prunet, N., & Jack, T. P. Flower Development in Arabidopsis: There Is More to It Than Learning Your ABCs. *Methods Mol. Biol.* **1110** 3-33 (2014).
2. Stewart, D., Graciet, E., & Wellmer, F. Molecular and regulatory mechanisms controlling floral organ development. *FEBS J.* (2016).
3. Heisler, M. G. *et al.* Patterns of auxin transport and gene expression during primordium development revealed by live imaging of the Arabidopsis inflorescence meristem. *Curr. Biol.* **15** (21), 1899-1911 (2005).
4. Reddy, G. V., & Meyerowitz, E. M. Stem-cell homeostasis and growth dynamics can be uncoupled in the Arabidopsis shoot apex. *Science*. **310** (5748), 663-667 (2005).
5. Yadav, R. K. *et al.* WUSCHEL protein movement mediates stem cell homeostasis in the Arabidopsis shoot apex. *Genes & Dev.* **25** (19), 2025-2030 (2011).
6. Reddy, G. V., Heisler, M. G., Ehrhardt, D. W., & Meyerowitz, E. M. Real-time lineage analysis reveals oriented cell divisions associated with morphogenesis at the shoot apex of Arabidopsis thaliana. *Development*. **131** (17), 4225-4237 (2004).
7. Chandler, J. W., Jacobs, B., Cole, M., Comelli, P., & Werr, W. DORNROSCHEN-LIKE expression marks Arabidopsis floral organ founder cells and precedes auxin response maxima. *Plant Mol. Biol.* **76** (1-2), 171-185 (2011).
8. Urbanus, S. L. *et al.* In planta localisation patterns of MADS domain proteins during floral development in Arabidopsis thaliana. *BMC Plant Biol.* **9** 5 (2009).
9. Hong, L. *et al.* Variable Cell Growth Yields Reproducible OrganDevelopment through Spatiotemporal Averaging. *Dev. Cell.* **38** (1), 15-32 (2016).

10. Roeder, A. H. *et al.* Variability in the control of cell division underlies sepal epidermal patterning in *Arabidopsis thaliana*. *PLoS Biol.* **8** (5), e1000367 (2010).
11. Heisler, M. G., & Ohno, C. Live-imaging of the *Arabidopsis* inflorescence meristem. *Methods Mol. Biol.* **1110** 431-440 (2014).
12. Tobin, C. J., & Meyerowitz, E. M. Real-Time Lineage Analysis to Study Cell Division Orientation in the *Arabidopsis* Shoot Meristem. *Methods Mol. Biol.* **1370** 147-167 (2016).
13. Smyth, D. R., Bowman, J. L., & Meyerowitz, E. M. Early flower development in *Arabidopsis*. *Plant Cell.* **2** (8), 755-767 (1990).
14. Prunet, N., Jack, T. P., & Meyerowitz, E. M. Live confocal imaging of *Arabidopsis* flower buds. *Dev. Biol.* (2016).
15. Fernandez, R. *et al.* Imaging plant growth in 4D: robust tissue reconstruction and lineaging at cell resolution. *Nat. Methods.* (2010).
16. Cutler, S. R., Ehrhardt, D. W., Griffiths, J. S., & Somerville, C. R. Random GFP::cDNA fusions enable visualization of subcellular structures in cells of *Arabidopsis* at a high frequency. *Proc. Natl. Acad. Sci. U S A.* **97** (7), 3718-3723 (2000).
17. Irish, V. F., & Sussex, I. M. Function of the *apetala-1* gene during *Arabidopsis* floral development. *Plant Cell.* **2** (8), 741-753 (1990).
18. Vernoux, T. *et al.* The auxin signalling network translates dynamic input into robust patterning at the shoot apex. *Mol. Systems Biol.* **7** 508 (2011).
19. Zhou, Y. *et al.* Control of plant stem cell function by conserved interacting transcriptional regulators. *Nature.* **517** (7534), 377-380 (2015).
20. Barbier de Reuille, P. *et al.* MorphoGraphX: A platform for quantifying morphogenesis in 4D. *Elife.* **4** 05864 (2015).
21. Barbier de Reuille, P. B., Robinson, S., & Smith, R. S. Quantifying cell shape and gene expression in the shoot apical meristem using MorphoGraphX. *Methods Mol. Biol.* **1080** 121-134 (2014).